

All you see is *all you need*

## ProActive® Implant System Box

- ✓ Few instruments
- ✓ Easy to clean and sterilize
- ✓ Surgical and prosthetic instruments in one box

### Features & benefits

**The Box** – Compact box where the lid is openable and removable with just one hand.

**The Trays** – Grommet-free modular platforms that keeps your instruments in place.


The grey surgical tray has the drill protocol, with NP and SP marking and room for extra surgical instruments for storage purpose.

The red prosthetic tray has the instruments arranged for ease of use for organization and storage.

**Depth Gauge Probe** – It has two sides with drill markings and different functions. One side is fitting in the preparation osteotomy made by the Ø 2.2 mm drill and has, apart from the depth markings, a ball end for inspection of the bone surface with tactile feel. The other side is for the preparation osteotomy made by the Ø 3.0 mm drill and larger to verify drill depth. The handle is shaped to improve ergonomic handling for the user.



Depth Gauge Probe

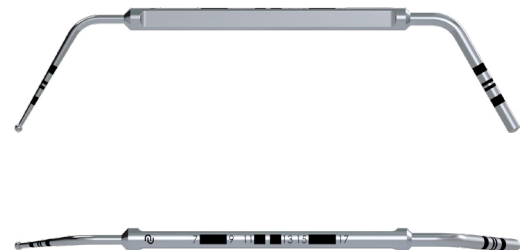


### Assortment

Art. No.	Description
51186	Neoss ProActive Implant System Box
51187	Neoss Surgical and Prosthetic Tray – ProActive Straight
51188	Neoss Surgical and Prosthetic Tray – ProActive Tapered
51189	Neoss Surgical and Prosthetic Tray – ProActive Edge
51190	Depth Gauge Probe


*All products might not be available in all markets.*

# Overview Surgical & Prosthetic Set-up

## ProActive® Edge Ø3.5 – 5.0 mm

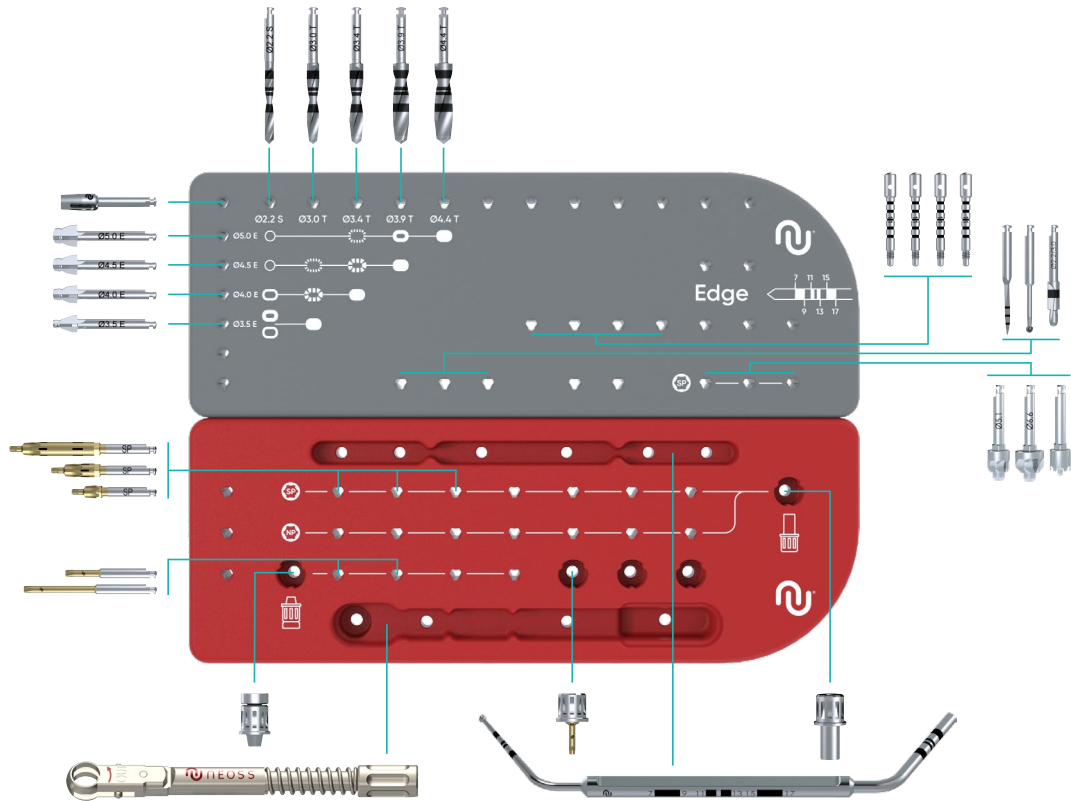

## ProActive® Straight Ø3.5 – 5.5 mm + ProActive® Straight, NP Ø3.25

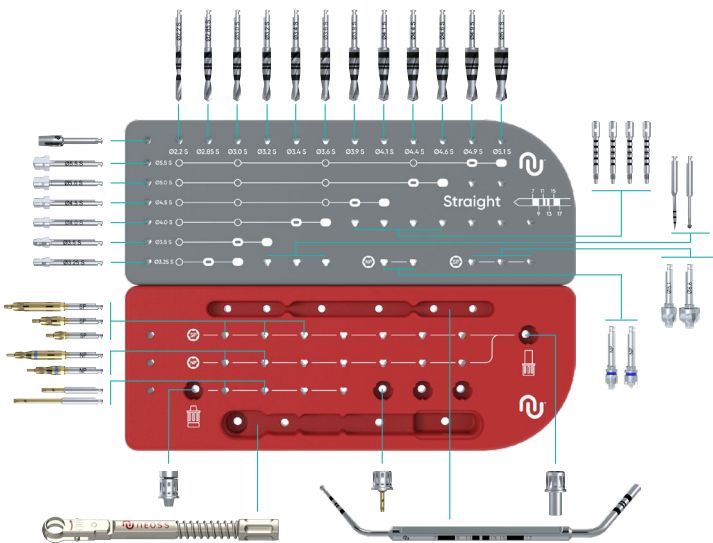

## ProActive® Tapered Ø3.5 – 5.5 mm

